



ERGONOMIC PIPETTING FOR PROBLEM LIQUIDS



| CONTAMINATION-FREE PIPETTING

Minimize the potential for pipette, sample, and human contamination with MICROMAN[®] E. Dispose of the capillary piston (CP) tip safely with a push of the button.

| RELIABLE RESULTS EVEN WITH PROBLEM LIQUIDS




Whether a viscous, volatile, hot, or cold sample, MICROMAN E can handle the job with accuracy and precision.

Have complete confidence in your results with the new volume-locking switch that reduces the risk of accidental volume adjustment.

| PROTECT YOUR PIPETTES

Positive-displacement technology makes it the ideal pipette for volatile liquids that can contaminate and corrode mechanical pipette pistons.

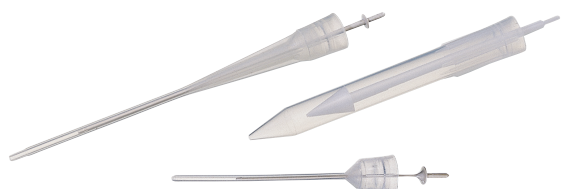

100% PROTECTION AGAINST CONTAMINATION

The disposable piston is located inside the capillary tip so users can be assured that no harmful samples will contaminate the pipette.

- Dye-free — no risk of chemical contamination.
- Made of pure polypropylene.
- Lightweight and slim CPs fit in even the narrowest of tubes.

MICROMAN® E			
MODEL	PART NUMBER	CAPILLARY PISTON	VOLUME (mL)
M10E	FD10001	CP10 CP10ST	1 5 10
M25E	FD10002	CP25 CP25ST	3 10 25
M50E	FD10003	CP50 CP50ST	20 50
M100E	FD10004	CP100 CP100ST	10 50 100
M250E	FD10005	CP250 CP250ST	50 100 250
M1000E	FD10006	CP1000 CP1000ST	100 500 1000

CAPILLARY PISTONS			
CAPILLARY PISTON	MICROMAN MODEL	PART NUMBER	QTY/BOX
CP10	M10E	F148412	192
		F148312	960
CP25	M25	F148112	200
		F148012	960
CP50	M50E	F148113	200
		F148013	960
CP100	M100E	F148414	192
		F148314	960
CP250	M250E	F148114	200
		F148014	960
CP1000	M1000E	F148560	182

STERILIZED CAPILLARY PISTONS			
CP10ST	M10E	F148413	192
		F148313	960
CP25ST	M25E	F148712	576
CP50ST	M50E	F148713	576
CP100ST	M100E	F148415	192
		F148315	960
CP250ST	M250E	F148714	576
CP1000ST	M1000E	F148180	182