



pipetman **NEO**[®]

A multichannel unlike any other!



4 MODELS OF PIPETMAN NEO MULTI WITH EXCLUSIVE FEATURES

| Model | Ordering reference | Volume Range |
|----------|--------------------|-----------------|
| 8 x 20 | F14401 | 2 µL to 20 µL |
| 12 x 20 | F14402 | 2 µL to 20 µL |
| 8 x 200 | F14403 | 20 µL to 200 µL |
| 12 x 200 | F14404 | 20 µL to 200 µL |

PATENTED PISTON DRIVE MECHANISM


EFFORT LESS TIP EJECTION

VERSATILE TIP-HOLDER



1 Innovative patented piston drive mechanism

Absolute consistency in volume distributed!

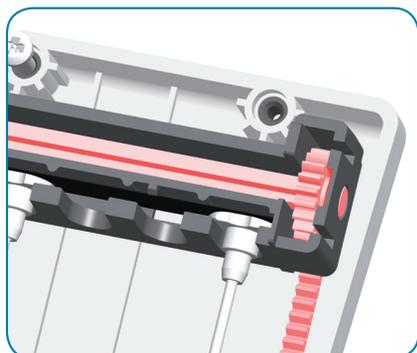


The new rack rail mechanism guarantees identical performance on each channel.

Thanks to an exclusive patented piston drive mechanism, the **pressure is applied consistently and simultaneously on each piston**, resolving the common downfall of multichannel pipettes: disparities between channels of the volume distributed.

2 Ingenious tip ejector technology

Best protection from Repetitive Strain Injuries (RSI)!

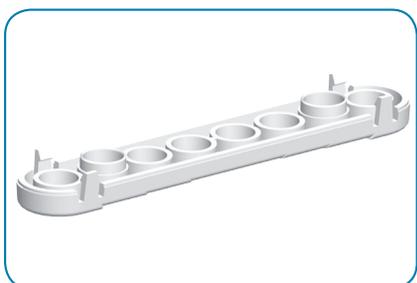


For more comfort, the ejector spacer has been designed to apply ejection forces over each tip sequentially.

This symetric system **drastically reduces tip ejection forces and risk of RSI**.

3 Unique tip-holder design

No need to purchase dedicated tip brand!

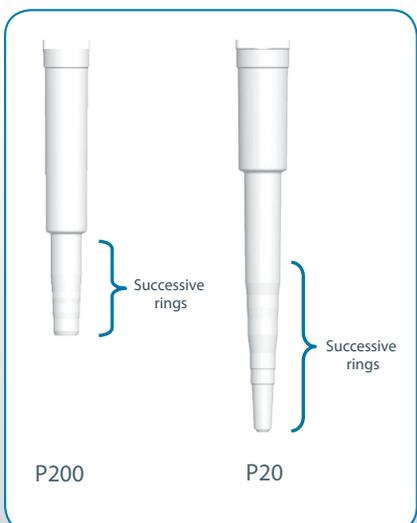

Successive rings have been positioned depending on the pipette volume range to guaranty the snug fit of tips throughout the protocols.

PIPETMAN DIAMOND® Tips are the perfect match to meet ISO's recommendations, but PIPETMAN Neo MULTI is also **fully compatible with the most commonly used tips on the market** such as: Axygen, Greiner, Starlab, Sarsstedt, etc.



www.gilson.com | sales@gilson.com | service@gilson.com | training@gilson.com

Gilson, Inc. World Headquarters | 3000 Parmenter Street | P.O. Box 620027
Middleton, WI 53562-0027 USA
Tel: (1) 800-445-7661 or (1) 608-836-1551 | Fax: (1) 608-831-4451

Gilson S.A.S. | 19, avenue des Entrepreneurs | BP 145, F-95400 Villiers-le-Bel, FRANCE
Tel: (33-1) 34 29 50 00 | Fax: (33-1) 34 29 50 20

 **GILSON®**

ISO9001 Certified

LT800504 © 06/2010 Gilson SAS. All rights reserved.
All trademarks are the property of Gilson Inc.
and its subsidiaries.